

crystalAligner v1.0 - Documentation


1 ABOUT

crystalAligner is a computer program for alignment of crystals in a scanning electron microscope. Given one or two crystal orientations (obtained from i.e. EBSD), one or two crystal directions and associated axes of the microscope coordinate system as alignment objectives and the rotational axes and limits of the microscope stage, the program optimizes the alignment of the crystal(s) with the microscope coordinate system by global optimization. The strength of *crystalAligner* is that it can utilize crystal symmetries and works with a heuristic global optimization algorithm. Even in the case of a standard microscope stage the alignment can therefore often be optimized to a satisfactory degree. A scientific treatment of the crystallographic aspects and the optimization procedure can be found in the associated publication [in Preparation].

2 PLATFORM AND INSTALLATION

The program was developed in MATLAB and requires the installation of the crystallographic texture analysis toolbox MTEX as well as the MATLAB global optimization toolboxes. The [global optimization toolbox](#) often comes preinstalled with MATLAB or can be obtained separately. If you are not having the toolbox installed you should receive an error message related to that. [MTEX](#) can be downloaded free of charge. The *crystalAligner* program has been tested with MTEX versions 4.2.1 to 5.1.1 and MATLAB version R2016b, R2015a and R2014a.

From version R2016b forward MATLAB changed the algorithm for linear constraints in *ga* (details [here](#)). Older versions of *ga* will experience problems of running into local minima (see figure below). While no ideal fix for this is possible, increasing the population size will mitigate this problem (try increasing the initial population size from 100 to 200 in examples E1 and E4). Multiobjective optimization problems using *gamultiobj* are not affected by this.

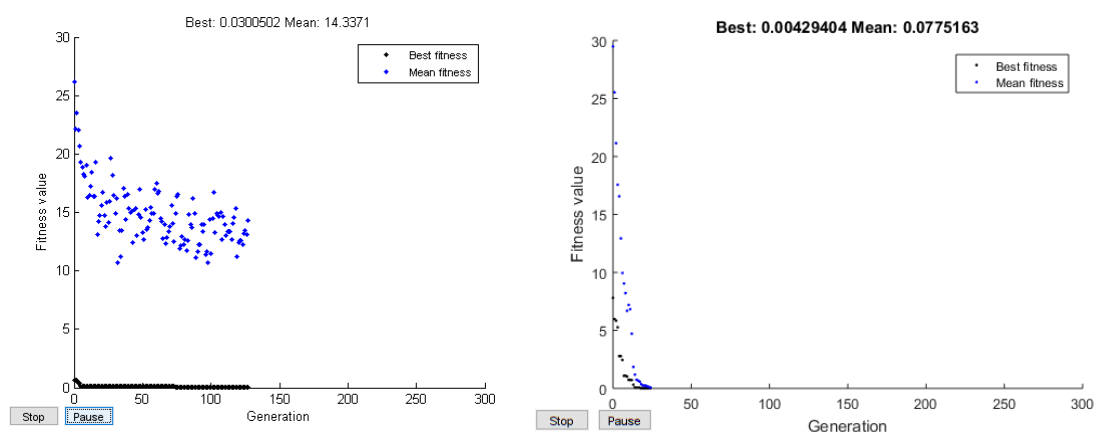

Comparison of the convergence of the mean and best fitness value in the MATLAB R2014a vs R2016b solution of problem 2 in example E4.

On an installed MATLAB system including the global optimization toolbox, MTEX can be installed by following these [instructions](#). You will have to navigate to your MTEX folder, right-click on it and choose 'Add to Path' and 'Selected Folder and Subfolders' to add the MTEX directory to your

MATLAB path. If your *crystalAligner.m* file is not in your MTEX folder, you will also have to add the path of your *crystalAligner.m* file to the MATLAB path in the same way. The program can be executed by entering 'crystalAligner' into the command window and pushing enter. It is however recommended to open the *crystalAligner.m* file in the editor, as the settings of the optimization problem should be checked and if necessary adjusted. From here the program can be executed by clicking on the play button.

3 DIFFERENT PARAMETER SETTINGS

All necessary input to the computer program is given in the section 'Initialization' in the top of the program. The meaning of each parameter is well documented in the program. A more detailed description of parameters is also treated in the associated publication [in Preparation]. In the following, some example parameter sets are presented. Keep in mind that the optimization algorithm is a genetic algorithm, meaning that the presented solutions are not going to be exactly identical for each execution of the code, but very similar (... practically speaking identical).

E1: ALIGNMENT OF A EITHER OF TWO EQUIVALENT CRYSTAL DIRECTIONS IN A CUBIC CRYSTAL

Two equivalent crystal direction in a cubic crystal with Euler orientation [233 145 33] are aligned in two separate single-objective optimization problems. The crystal represents a ferrite grain in a steel sample. The optimization objective is to align the z-axis of the microscope with either of the $\langle 111 \rangle$ or $\langle 100 \rangle$ equivalent directions. The stage axes represent a standard stage of an FEI scanning electron microscope with chosen soft-limits in the tilt axis of $0 - 20^\circ$. The genetic algorithm is set up to have a population size of 100 individuals. Optional text output for estimating the trench sizes for a FIB lift-out at an inclined surface are deactivated. It is found that the $\langle 100 \rangle$ can be perfectly aligned while the optimal $\langle 111 \rangle$ alignment still deviates about 20° from ideal alignment, caused by the significant constraints by the SEM stage rotation limits.

Settings for running this example:

```
% *** General settings ***
%Crystal
crys.o = [313 15 137];
crys.cs = {'cubic'};
crys.ss = 'orthorhombic';
%Alignment Objective(s)
axs.align = [zvector];
dir.ax{1} = [1 1 1; 1 0 0];
axs.sym = [ 1 ];
%SEM stage definition
axs.rot = [xvector; zvector];
optim.LB = [ 0 -180 ];
optim.UB = [ 20 180 ];
axs.sign = [ 1 -1 ];
axs.order = [ 2 1 ];

% *** Genetic algorithm - optimization settings
%Genetic algorithm
optim.popSz = 100;
optim.funcTol = 0.01;
optim.maxStallGen = 10;
optim.iterOut = 0;
%Multiobjective genetic algorithm settings
```

```

optim.wghtFac = [1,1];
optim.multiCore = 0;
optim.hybridFcn = 0;
optim.autoSol = 1;

% *** Optional - FIB liftout calculations
FIB.mode = 0;
FIB.trench.ang = 52;
FIB.trench.z = 15;
FIB.axs.tilt = 1;
FIB.axs.rot = 2;
% *** Output
optim.plot = 1;

```

Expected result output in command window:

```

*****
Optimization problem 1 - Parallel alignment of:
  -> m-3m 111 direction with microscope Z-axis
*****
Optimization terminated: average change in the fitness value less than
options.FunctionTolerance.
-----
*** Optimization results ***
Optimal parallel alignment of microscope Z-axis with crystal direction 111:
  -> Rotation around microscope X-axis: 19.9°
  -> Rotation around microscope Z-axis: -39.3°
  -> Deviation from ideal alignment in Z-axis: 20.0 °
-----

*****
Optimization problem 2 - Parallel alignment of:
  -> m-3m 100 direction with microscope Z-axis
*****
Optimization terminated: average change in the fitness value less than
options.FunctionTolerance.
-----
*** Optimization results ***
Optimal parallel alignment of microscope Z-axis with crystal direction 100:
  -> Rotation around microscope X-axis: 15.0°
  -> Rotation around microscope Z-axis: 133.0°
  -> Deviation from ideal alignment in Z-axis: 0.0 °
-----

```

E2: ALIGNMENT OF TWO EQUIVALENT CRYSTAL DIRECTIONS IN A CUBIC CRYSTAL

This case is equivalent to the application example 1 in the publication [in Preparation]. Here, two equivalent crystal directions in a cubic crystal with Euler orientation [61 42 9] are aligned. The crystal represents a β -phase grain in a near- β Ti alloy. The optimization objectives are to align the y-axis of the microscope with any of the $\langle 113 \rangle$ equivalent directions and the z-axis of the microscope with any of the $\langle 110 \rangle$ equivalent directions. The stage axes represent a standard stage of an FEI scanning electron microscope with chosen soft-limits in the tilt axis of $0 - 20^\circ$. The genetic algorithm is set up to have a population size of 100 individuals and no weighting factors applied to the objectives for the choice of the optimal solution. Optional text output for estimating the trench sizes for a FIB lift-out at an inclined surface and plotting of graphs are activated.

Settings for running this example:

```

% *** General settings ***
%Crystal
crys.o = [61 42 9];
crys.cs = {'cubic'};
crys.ss = 'orthorhombic';
%Alignment Objective(s)
axs.align = [yvector; zvector];
dir.ax{1} = [1 1 3];
dir.ax{2} = [1 1 0];
axs.sym = [ 1 1 ];
%SEM stage definition
axs.rot = [xvector; zvector];
optim.LB = [ 0 -180 ];
optim.UB = [ 20 180 ];
axs.sign = [ 1 -1 ];
axs.order = [ 2 1 ];

% *** Genetic algorithm - optimization settings
%Genetic algorithm
optim.popSz = 100;
optim.funcTol = 0.01;
optim.maxStallGen = 10;
optim.iterOut = 0;
%Multiobjective genetic algorithm settings
optim.wghtFac = [1,1];
optim.multiCore = 0;
optim.hybridFcn = 0;
optim.autoSol = 1;

% *** Optional - FIB liftout calculations
FIB.mode = 1;
FIB.trench.ang = 52;
FIB.trench.z = 15;
FIB.axs.tilt = 1;
FIB.axs.rot = 2;

% *** Output
optim.plot = 1;

```

Expected result output in command window:

```

*****
Optimization problem 1 - Parallel alignment of:
  -> m-3m 113 direction with microscope Y-axis
  -> m-3m 110 direction with microscope Z-axis
*****

```

Optimization terminated: average change in the spread of Pareto solutions less than options.FunctionTolerance.

```

*****
Optimal solution found by shortest distance to origin
under consideration of weighting factors 1 and 1 for alignment with
microscope Y-axis and Z-axis

```

```

-----
*** Optimization results ***
Optimal parallel alignment of microscope Y-axis with m-3m crystal direction
113 and microscope Z-axis with m-3m crystal direction 110:
  -> Rotation around microscope X-axis: 6.9°
  -> Rotation around microscope Z-axis: 132.3°

```

```

-> Deviation from ideal alignment in Y-axis: 0.1 °
-> Deviation from ideal alignment in Z-axis: 0.5 °
-----
*** FIB - FEI Helios ***
Apply:
-> Relative rotation of 132.3°
-> Tilt at lift-out position: 6.8°
-> Tilt at trenching position: 58.8°
-> Alternative: Tilt at trenching position: 45.2° + 180° relative
rotation
Trench lengths for 15.0 µm trench depth (z) and 52.0° trench angle:
-> Trench length (y) at 'downhill position': 16.6 µm
-> Trench length (y) at 'uphill position': 22.7 µm
-----

```

E3: ALIGNMENT OF TWO CRYSTAL DIRECTIONS IN TWO CRYSTALS OF DIFFERENT STRUCTURE AND ORIENTATION

This case is equivalent to the application example 2 in the publication [in Preparation]. Here, crystal directions in a cubic and an orthorhombic crystal with Euler orientations [261 43 28] and [175 20 102] are aligned. They represent a β -parent grain and an α'' martensitic variant in a near- β Ti deformation microstructure. The optimization objectives are to align the x-axis of the microscope with any of the $\langle 110 \rangle$ equivalent directions in β and the z-axis of the microscope with the [001] axis of the martensitic α'' . The stage axes represent a standard stage of an FEI scanning electron microscope with chosen soft-limits in the tilt axis of 0 – 20°. The genetic algorithm is set up to have a population size of 500 individuals and no weighting factors applied to the objectives for the choice of the optimal solution. Optional text output for estimating the trench sizes for a FIB lift-out at an inclined surface is deactivated while plotting of graphs is activated.

Settings for running this example:

```

% *** General settings ***
%Crystal
crys.o = [261 43 28; 175 20 102];
crys.cs = {'cubic', 'orthorhombic'};
crys.ss = 'orthorhombic';
%Alignment Objective(s)
axs.align = [xvector; zvector];
dir.ax{1} = [0 1 1];
dir.ax{2} = [0 0 1];
axs.sym = [ 1 0 ];
%SEM stage definition
axs.rot = [xvector; zvector];
optim.LB = [ 0 -180 ];
optim.UB = [ 20 180 ];
axs.sign = [ 1 -1 ];
axs.order = [ 2 1 ];

% *** Genetic algorithm - optimization settings
%Genetic algorithm
optim.popSz = 500;
optim.funcTol = 0.1;
optim.maxStallGen = 10;
optim.iterOut = 0;
%Multiobjective genetic algorithm settings
optim.wghtFac = [1,1];
optim.multiCore = 0;

```

```

optim.hybridFcn = 0;
optim.autoSol = 1;

% *** Optional - FIB liftout calculations
FIB.mode = 0;
FIB.trench.ang = 52;
FIB.trench.z = 15;
FIB.axs.tilt = 1;
FIB.axs.rot = 2;

% *** Output
optim.plot = 1;

```

Expected result output in command window:

```

*****

Optimization problem 1 - Parallel alignment of:
  -> m-3m 011 direction with microscope X-axis
  -> mmm 001 direction with microscope Z-axis
*****

Optimization terminated: average change in the spread of Pareto solutions
less than options.FunctionTolerance.

*****

Optimal solution found by shortest distance to origin
under consideration of weighting factors 1 and 1 for alignment with
microscope X-axis and Z-axis

-----
*** Optimization results ***
Optimal parallel alignment of microscope X-axis with m-3m crystal direction
011 and microscope Z-axis with m-3m crystal direction 001:
  -> Rotation around microscope X-axis: 19.5°
  -> Rotation around microscope Z-axis: 6.4°
  -> Deviation from ideal alignment in X-axis: 6.6 °
  -> Deviation from ideal alignment in Z-axis: 3.9 °
-----

```

E4: EXTENSION TO EXAMPLE E1 BY ADDING A SUBSTAGE

In example E1 it was found that the $\langle 100 \rangle$ can be perfectly aligned while the optimal $\langle 111 \rangle$ alignment still deviates about 20° from ideal alignment, caused by the significant constraints by the SEM stage rotation limits. Let's investigate the effect of adding a sub-stage on the achievable alignment for $\langle 111 \rangle$. Compared to the setup for example E1 a rotation axis around the microscope y-axis was added (axs.rot). The rotation limits for this axis are defined as $-45 - 45^\circ$. Here it is important to consider the hierarchy of the rotational axes: The sub-stage represents the 1st tilt rotation which should be aligned with the microscope y-axis as an initial position, the microscope z-axis the 2nd rotation and the microscope x-axis the 3rd rotation (axs.order = [3 1 2]). The additional degree of freedom compared to the setup in example E1 leads to perfect alignment of the $\langle 111 \rangle$ direction with the microscope z coordinate axis and is reflected in that repetitive execution of the code will produce different sets of stage-rotations leading to the same angular deviation. As an additional example one could try to set the limits of the substage to i.e. $45 - 45^\circ$, which would simulate the use of a fixed 45° pre-tilt sample holder.

Settings for running this example:

```

% *** General settings ***
%Crystal
crys.o = [313 15 137];
crys.cs = {'cubic'};
crys.ss = 'orthorhombic';
%Alignment Objective(s)
axs.align = [zvector];
dir.ax{1} = [1 1 1; 1 0 0];
axs.sym = [ 1 ];
%SEM stage definition
axs.rot = [xvector; yvector; zvector];
optim.LB = [ 0 -45 -180 ];
optim.UB = [ 20 45 180 ];
axs.sign = [ 1 1 -1 ];
axs.order = [ 3 1 2 ];

% *** Genetic algorithm - optimization settings
%Genetic algorithm
optim.popSz = 100;
optim.funcTol = 0.01;
optim.maxStallGen = 10;
optim.iterOut = 0;
%Multiobjective genetic algorithm settings
optim.wghtFac = [1,1];
optim.multiCore = 0;
optim.hybridFcn = 0;
optim.autoSol = 1;

% *** Optional - FIB liftout calculations
FIB.mode = 0;
FIB.trench.ang = 52;
FIB.trench.z = 15;
FIB.axs.tilt = 1;
FIB.axs.rot = 2;
% *** Output
optim.plot = 1;

```

Expected result output in command window:

```

*****
Optimization problem 1 - Parallel alignment of:
-> m-3m 111 direction with microscope Z-axis
*****

```

Optimization terminated: average change in the fitness value less than options.FunctionTolerance.

```

-----
*** Optimization results ***
Optimal parallel alignment of microscope Z-axis with crystal direction 111:
-> Rotation around microscope X-axis: 14.5°
-> Rotation around microscope Y-axis: 37.4°
-> Rotation around microscope Z-axis: -116.3°
-> Deviation from ideal alignment in Z-axis: 0.0 °
-----

```

```

*****
Optimization problem 2 - Parallel alignment of:
-> m-3m 100 direction with microscope Z-axis
*****

```

Optimization terminated: average change in the fitness value less than options.FunctionTolerance.

```
-----  
*** Optimization results ***  
Optimal parallel alignment of microscope Z-axis with crystal direction 100:  
-> Rotation around microscope X-axis: 2.2°  
-> Rotation around microscope Y-axis: 14.8°  
-> Rotation around microscope Z-axis: 51.2°  
-> Deviation from ideal alignment in Z-axis: 0.0 °  
-----
```

4 TROUBLESHOOTING

4.1 ERROR 'NUMBER OF AXIS-CRYSTALDIRECTION PAIRS OF ALIGNMENT OBJECTIVES NOT EQUAL'

You have probably given more microscope axes for alignment (axs.align) than number of crystal directions for alignment in (dir.ax) or vice versa. The number of alignment axes, crystal directions and symmetry flags in the initialization part 'Alignment Objective(s)' have to have equal dimensions as demonstrated in the examples.

4.2 ERROR 'NUMBER OF CRYSTALSYSTEM-ORIENTATION PAIRS NOT EQUAL'

In the initialization section 'General settings -> Crystal' it is important to list an equal amount of crystal structures (crys.cs) and crystal orientations (crys.o).

4.3 ERROR WHEN STARTING AND INITIALIZING MTEX

The script attempts to automatically start and initialize MTEX. If an error occurs try to restart MATLAB and start and initialize MTEX manually by executing 'startup_mtex' before running the *crystalAligner* program.

4.4 OUTPUT 'LINEAR CONSTRAINTS ARE NOT SATISFIED WITHIN CONSTRAINT TOLERANCE'

This is an error that occurs for single-objective optimization (examples E1 and E4) and is related to an implementation of the *ga* algorithm in MATLAB versions older than R2016b (see textbox on page 1). From R2016b on the *ga* algorithm will strictly fulfil the set linear constraints, even if the achieved alignment is highly non-optimal. It likely means that you ran into a local minimum; while installing a newer version of MATLAB would be the best way to go a quick-and-dirty way of avoiding this is to increase the population size. This will slow down the code performance.

4.5 OTHER ISSUES

If there are other issues, comments or suggestions you would like to share with me you can contact me via contactnospam@fniessen.com (remove the nospam to make this email address work) and I will be happy to help.